

Technical Information

Primax Roll Fix

X-ray Film Fixer for Roller Machines (liquid)

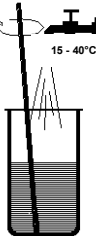
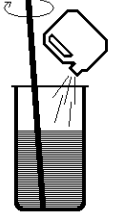
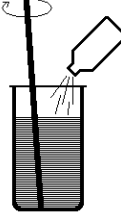
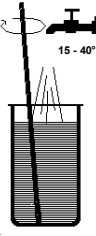
Description:

- X-ray Film fixer for Automatic Processing in processors of roller transport type with fixing times at 15 s or more
- Because of its excellent hardening properties and fast working fixer components this product is suitable for use in fast going machines.
- application in conjunction with Developer Primax Roll DEV
- all liquid concentrates, two parts, for medium to high production volume

Replenishment :

For fixing of 1 m² X-ray film **600 to 800 ml** replenisher are necessary at standard conditions . Generally replenishment occurs automatically by adding solution depending on film amount.

Mix Instructions :

Water	Part A	Part B	Water
			
10,0 l	5,0 l	1,25 l	3,75 l

First dissolve part A completely in water, afterwards add part B under permanent stirring. Finally add water to get 20 l ready-to-use product.

Assortment :

Order number	Product Name	Mix Size for	Contents per unit (concentrates)	Capacity for ... m ² film
209505	Primax Roll FIX	2 x 20 l	Part A 2 x 5,0 l Part B 2 x 1,25 l	50 - 65

Storage and Handling :

- Until the preparation the concentrates has to be kept in the closed original package and stored at temperatures of 10 to 25 °C.
- X-ray Film chemicals contain components which may cause irritations in case of contact with the skin.
- In case of contact with the eyes it has to be rinsed thoroughly with plenty of running cold water by opened eyelid and a medical doctor has to be consulted.